# Jordan CLEAN TRAVEL

## Electrical Toothbrush

Bruksanvisning

Brugsanvisning

Käyttöohje

Instruction manual

TBPL-100GY

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[www.wilfa.fi](http://www.wilfa.fi)

## JORDAN CLEAN ELECTRIC TOOTHBRUSH

The Jordan Clean series is a new generation of electric toothbrushes that through SoftClean™ Technology delivers the right combination of oral hygiene features. Inspired by a Scandinavian design approach for simplicity and user-friendliness, the series gives an effective clean from an electric toothbrush and the ease of use of a manual one.

### SAFETY

Please note the following instructions when using the appliance:

- This appliance is intended for domestic use only.
- This appliance is designed solely for cleaning and caring for human teeth. All other uses are expressly prohibited.
- Do not use if the appliance or attachments is visibly damaged.
- When the appliance is being charged, on no account should this or the charging unit come into contact with water or any other liquid.
- Keep the appliance out of direct sunlight, excessive heat and open flames.
- Do not drop the appliance.
- Keep the appliance out of the reach of children.
- Do not open the appliance under any circumstances.
- No warranty claims are accepted for damage caused by improper handling.
- Cleaning and user maintenance shall not be made by children without supervision.
- This appliance can be used by children aged from 8 years and above. And persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- The Electric Toothbrush can only be charged by the charging unit provided with the appliance.

### SPECIAL FEATURES:

- **Softclean™ technology** gives you a unique combination of highly effective bristles, and a smaller head for effective cleaning in a effective way. Each head is divided into three cleaning zones with specialized efficiencies per zone.
- **2 modes** with high oscillating speed at about 8800 rpm at hi speed for increasing cleaning performance and soft brushing at low speed of approximately 7200 rpm. Hold the switch for 1s for selecting the desired Hi/Lo Speed.
- **2-min timer** with click signal for each 30 sec brushing time for reminding you to move to different teeth area. It will click 2 times at 30/60/90 sec. and click 4 times at 2 minutes (recommended brushing time). The brush will stop automatically after 10 minutes for protecting your teeth and preserve battery capacity by accidental switch-on when travelling.
- **Low power indication** when battery capacity goes down to about 20%, the indicator blinks quickly to remind you to recharge the brush in order to keep full power.
- **Built-in high capacity Lithium battery for long use.** It can run 120min which gives 1 month based on 2 times brushing a day and 2 minutes brushing each time.


- **Over-charge and over-discharge Protection** circuit is equipped to preserve the battery life time. When fully charged the toothbrush will automatic cut off. When running low on battery the toothbrush will automaticly stop.

## SPECIFICATIONS

Maximum operation time is approximately 120 minutes/4 weeks when fully charged. When rechargeable batteries are running low, the indicator will blink to tell you to recharge the unit immediately. If you continue to use the unit for some time without recharging it, the unit may stop automatically. This is our battery protection feature to preserve the batteries life time. You need to recharge immediately.

Power supply:

Jordan USB cable

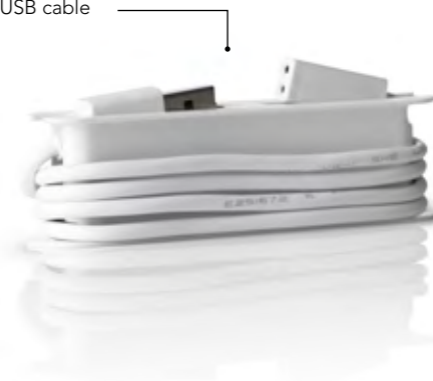MICRO USB adapter

Appliance:

Lithium-ion 3.7V 600mA  
charged by DC5V, 800mA.

## CHARGING THE UNIT

The toothbrush can be charged with the supplied cable that fits into a USB port. It is also possible to use any MICRO USB / USB cable with the supplied adapter.

**IMPORTANT:** Remember to close the lid in the bottom of the toothbrush after charging.

The indicator will blink with light/shade (breathing mode) alternately, when the power unit is insert to charging base which is plugged-in. It will light up WHITE when batteries are fully charged (18 hours). The power unit may warm up somewhat during recharging. This has no significance on the unit's operation. Over-charging is impossible.

## PRESERVING THE RECHARGEABLE BATTERIES

**To maximize the capacity and life of the rechargeable batteries, they should be fully charged before first use and then fully discharged through use. Then recharge to full capacity again. This discharging/charging procedure should be done every six months. The Electric Toothbrush also has built-in LOW BATTERY INDICATION feature. If the battery power often runs too low this will shorten the battery life. The indicator will blink to remind you to charge. If you don't charge another "Automatic Shut off system" will stop the brush. You need to fully recharge before using again.**

## BEFORE USE

Before using your Jordan Clean toothbrush for the first time, charge the unit full. After the power unit has been fully charged, attach the brush head to the power unit. The brush head must face the same direction as the switch after twist in. Jordan Clean has only one switch by pressing for the desired functions. In case you put the power unit into the charging base when it is working, it will automatically stop.

## REMOVAL OF RECHARGEABLE BATTERY

This toothbrush is using Lithium-ion rechargeable battery. When the toothbrush reaches the end of its service life, the battery must be removed from the appliance before it is discarded. The appliance must be disconnected from the mainpower supply before removing the battery. The battery must be disposed of safely. Use a flathead screwdriver to remove the internal body, then remove the battery and dispose of it via a return and collection system or contact the retailer from whom the product was purchased. They operate a take-back system for safe recycling.

## PROTECTING THE ENVIRONMENT

Correct disposal of this product.

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from

uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

**WILFA AS, Industriveien 25, 1481 Hagan, Norway**

[wilfa.com](http://wilfa.com)